




# FluoroMesh

## Fluoroscopy Mesh Phantom

A routine test object designed to provide an ongoing check of the resolution of fluoroscopy systems. Image quality is measured simply by counting the number of mesh patterns resolved in the image. An ongoing record of these numbers will reveal any trend towards deterioration in imaging performance.

Available in three resolution ranges:

- **FLUOROMESH 16-60**  
contains 8x mesh wedges with a range of 16 to 60 lines/inch
- **FLUOROMESH 30-100**  
contains 8x mesh wedges with a range of 30 to 100 lines/inch
- **FLUOROMESH 60-150**  
contains 8x mesh wedges with a range of 60 to 150 lines/inch

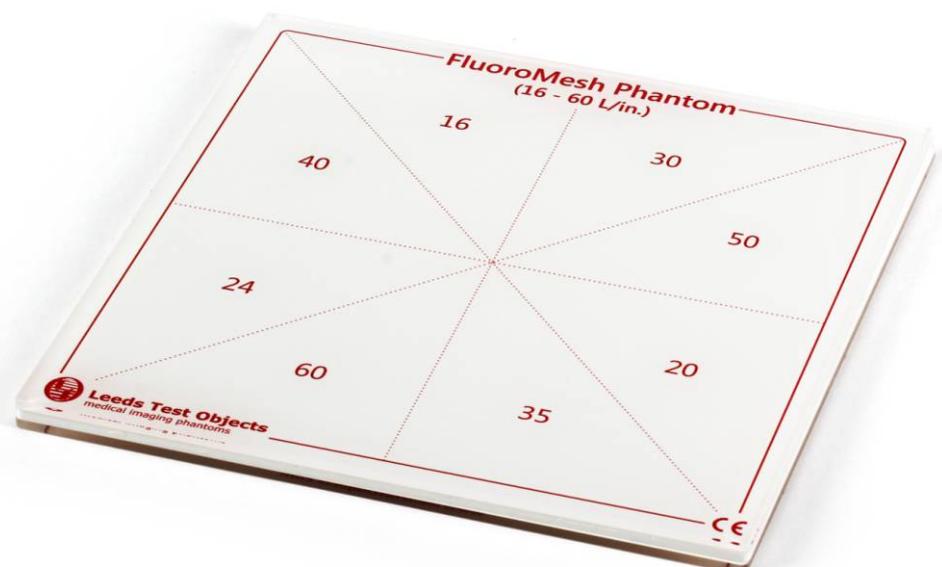



# FluoroMesh

Fluoroscopy

## Product X-rays

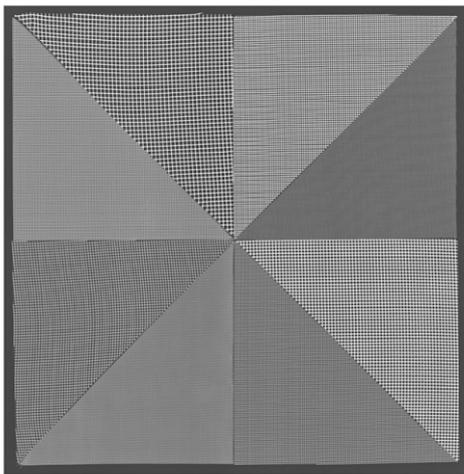


fig. 1 FluoroMesh 16-60 x-ray

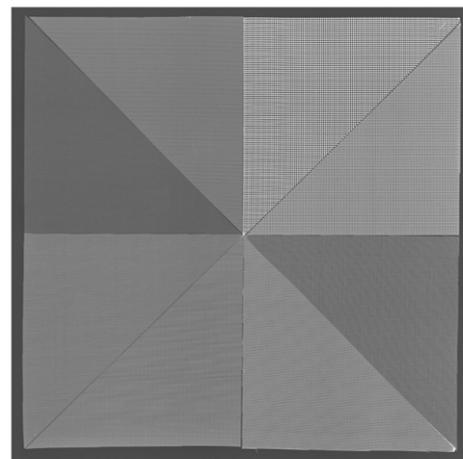


fig. 2 FluoroMesh 30-100 x-ray

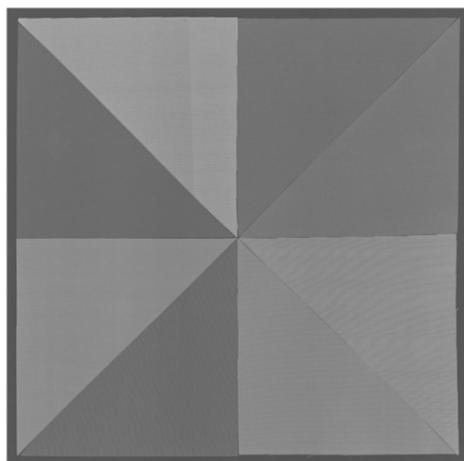


fig. 3 FluoroMesh 60-150 x-ray

[www.leadstestobjects.com](http://www.leadstestobjects.com)