



# TO PAN


## Dental Panoramic Phantom

A test object that allows the complex focal trough of the panoramic x-ray unit to be quantified, as well as the capacity for spatial resolution testing using line-pair test patterns (available separately).

This allows full assessment of parameters affecting image quality.  
Can be used in film and digital systems.  
Includes tripod for universal positioning.

### TO PAN Enables:

- Quantification of the focal trough and image layer
- Accurate determination of the bite offset
- Quantification of the symmetry of the scan
- Measurement of the incident beams angles
- Simple visualisation of faults
- Spatial resolution testing

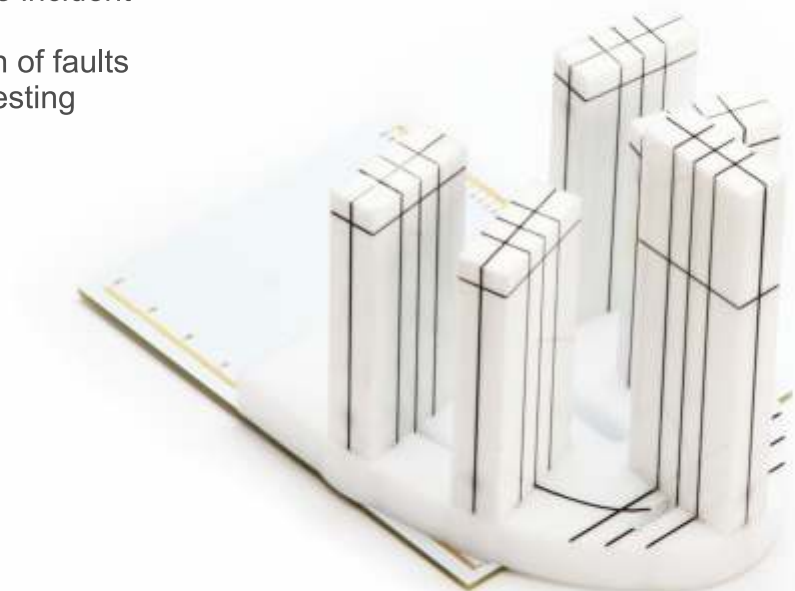

For more information visit us at [www.leadstestobjects.com](http://www.leadstestobjects.com)



# TO PAN

## Dental Panoramic Phantom



fig. 1 Radiograph of TO PAN



[www.leadstestobjects.com](http://www.leadstestobjects.com)